

# ELDERLY Safety Alert

## INCLUDES;

Main Unit (Receiver)	x1
Calling Units (Transmitter)	x2
Adapter 6V, 200mA ( <i>Optional</i> )	x1

## FEATURES

- A wireless signaling device with easy operation and functional design
- Great product to ensure the elderly get immediate attention at the push of a button
- Main unit can be powered either by 4pcs AA (UM-3) batteries (not included) or DC 6V adapter.
- 3 levels of alarm sound volume adjustable (High, Medium, Low)
- Low battery indicator
- 200 feet (open site) operation range
- Main Unit can be used in any room

## BATTERY INSTALLATION

### MAIN UNIT

1. Make sure the ON/OFF switch is in OFF position
2. Slide down the front cover from Main Unit
3. Follow the polarity markings and install 4pcs AA (UM-3) batteries (not included) into the battery compartment.


### CALLING UNIT

1. Unscrew the back cover and follow the polarity instructions to insert one fresh CR2032 (included) (see Drawing I)
2. Replace the cover and screw securely back on
3. After installing the battery, you will need to SET the Call Unit as instructed below



## HOW TO SET THE CALL UNIT?

1. After installing batteries, switch the Main Unit to ON. You will hear one beep and the RED LED will be ON to show the Main unit is properly functioning.
2. Slide down the Main Unit front cover and press the Learning-Code button inside CONTINUOUSLY (see Drawing II)

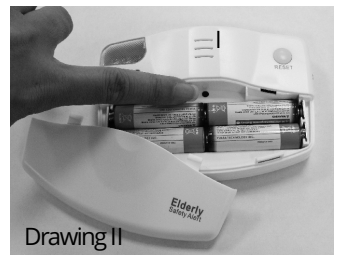

# ELDERLY Safety Alert

3. You will hear one beep and the Yellow/Green LED will illuminate.
4. Press the Learning-Code button and the Call unit simultaneously
5. You will hear one beep to indicate it is SET for Call unit 1
6. Repeat the procedures 2 ~ 5 to set the other Call units

## HOW TO OPERATE THE ELDERLY SAFETY ALERT

1. After SETTING the Call unit, press the Call unit button, the Main Unit will respond with a series of three-tone sounds and a flashing red LED
2. Press the Main unit RESET button to turn OFF the calling sound & RED LED
3. Different colored Yellow/Green LED will distinguish which Call unit (I or II) is sending the call message
4. Yellow LED is for Call unit 1
5. Green LED is for II Call unit 2

## POWER INDICATOR

### *Main Unit*

- If the battery power is low, the red LED will flash every 10-seconds
- Replace fresh 4pcs of AA (UM-3) batteries

### *Calling Unit*

- If the red LED becomes dim or the Calling Unit operating distance becomes shorter, the Calling Unit battery power is becoming low
- Replace the CR2032 battery
- After installing the Call unit battery, you have to SET the Call Unit again

## VOLUME ADJUSTMENT

- 3 alarm sound volume level (High, Medium, Low) adjustable
- SHORT press the Learning Code button and you will hear a corresponding alarm sound, SHORT press again to move to the next sound level
- Choose the suitable sound level and release the button to set the sound volume

## NOTE:

- The learning code has been memorized in the Main Unit, so there is no need to SET (learn code) after Main unit battery replacement. Only the Call unit battery replacement needs to be SET again (learn code again)
- The call sound will continue until the Main unit RESET button is pressed